



Microscopes for Flat Panel Display and Large Scale Integration Inspection

ECLIPSE L300N/L300ND L200N/L200ND

Microscopes for Flat Panel Display
and Large Scale Integration Inspection



ECLIPSE



L300N

For ø300 mm wafer/
Episcopic optical contrast



L300ND

For 17-inch FPD/
Episcopic and Diascopic optical contrast



L200N

For ø200 mm wafer/
Episcopic optical contrast



L200ND


For ø200 mm wafer/
Episcopic and Diascopic optical contrast

Enhanced observation performance and operation


Epi-fluorescence observation widens inspection range—including 365 nm UV excitation

*L300N/L300ND/L200ND only


- Highly beneficial when inspecting semiconductor resist residues and organic electroluminescence displays.
- Various observation methods such as brightfield, darkfield, simple polarizing, and DIC are possible on all models.
- With the L300ND/L200ND, diascopic illumination capability adds the illumination through transparent substrates.

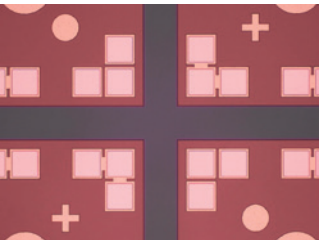

Brightfield observation of wafer pattern

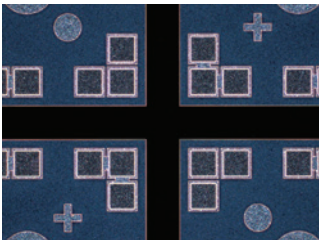

Darkfield observation

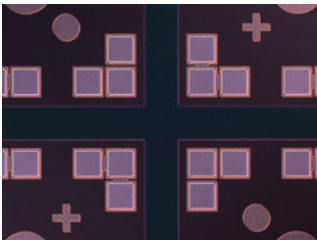


DIC observation



Epi-fluorescence observation of organic substance on wafer

Front operation with easy access


- Minimizes fatigue during lengthy observations, maintaining a safer operator distance from the sample

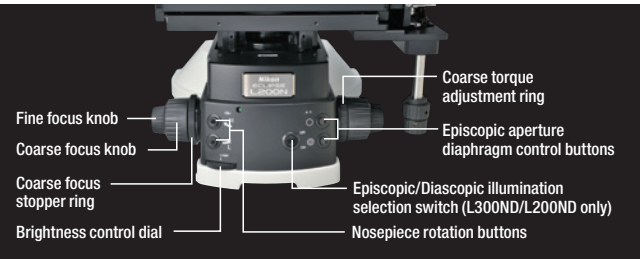

Fine focus knob
Coarse focus knob
Coarse focus stopper ring
Brightness control dial
Coarse torque adjustment ring
Episcopic aperture diaphragm control buttons
Episcopic/Diascopic illumination selection switch (L300ND/L200ND only)
Nosepiece rotation buttons

Target for easier focusing

- Insert a focusing target in the optical path to easily focus on low-contrast samples, such as bare wafers.

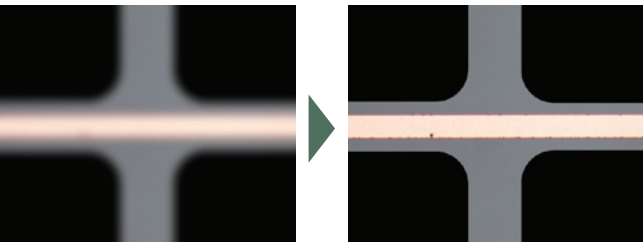


Stronger safeguard against contamination

- Antistatic coatings applied to the body, stage, eyepiece tube and other various controls
- Prevents damage to samples and contributes to higher yields

Observation at optimum eyepoint level

- Ultra-wide 25-mm field of view and eyepiece angle adjustment between 0 ° and 30 °
- Operators can adjust eyepoint level to ensure a comfortable viewing position
- L2-TT4A eyepiece tube enables centering and focus adjustment of the vertical tube.

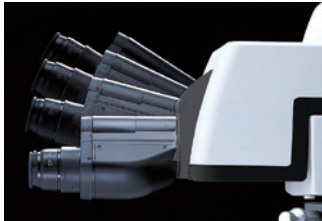

Fixed-position X-Y fine movement control

- Allows for stage movements and focusing to be carried out with ease
- An optional stage handle adapter is also available to enable finer stage control



The X-Y fine movement controls are positioned close to the operator.

Illumination

LED

Compact LED illuminators are power saving and achieve long life.



LV-LL LED Lamphouse

Intensilight

- Motorized mercury precentered fiber illuminator for epi-fluorescence observation, with variable light intensity and shutter control, provide excellent flexibility. Lamp centering and focus adjustment are not necessary.

Filter blocks

For epi-fluorescence observation

- EPI-FL UV-2A
- EPI-FL V-2A
- EPI-FL BV-2A



*Only one cube is attachable.

Accessories

Nikon's CFI60 optical systems are highly evaluated for their unique concept of high NA combined with a long working distance. These lenses have been developed further and evolved achieving the apex in long working distance specifications, correct chromatic aberration, and an optimized lens weight.

Objective lenses

Standard objective lenses

TU Plan Fluor Series

EPI/BD 5x/10x/20x/50x/100x

Enable brightfield, darkfield, simple polarizing, sensitive polarizing, differential interference, and epi-fluorescence observations with just one lens. Achieves superior chromatic aberration performance with long working distance for all magnifications to adapt to any application.



*Brightfield observation (EPI) objective lens

Model	Magnification	NA	Working Distance(mm)
TU Plan Fluor EPI (brightfield type)	5x	0.15	23.5
	10x	0.30	17.5
	20x	0.45	4.5
	50x	0.80	1.0
	100x	0.90	1.0
TU Plan Fluor BD (brightfield/ darkfield type)	5x	0.15	18.0
	10x	0.30	15.0
	20x	0.45	4.5
	50x	0.80	1.0
	100x	0.90	1.0

Long working distance objective lenses

TU Plan ELWD Series

EPI/BD 20x/50x/100x

With the phase Fresnel lenses, these objective lenses enable long working distances while offering higher level chromatic aberration correction than conventional objective lenses. This improves operability for samples with different heights.



*Brightfield observation (EPI) objective lens

Model	Magnification	NA	Working Distance(mm)
TU Plan EPI ELWD (brightfield type)	20x	0.4	19.0
	50x	0.6	11.0
	100x	0.8	4.5
TU Plan BD ELWD (brightfield/ darkfield type)	20x	0.4	19.0
	50x	0.6	11.0
	100x	0.8	4.5

Low-magnification objective lenses

T Plan EPI

EPI 1x/2.5x

Model	Magnification	NA	Working Distance(mm)
T Plan EPI (brightfield type)	1x	0.03	3.8
	2.5x	0.075	6.5

Apochromatic objective lenses

TU Plan Apo Series

EPI/BD 50x/100x/150x

By using phase Fresnel lenses, these objective lenses achieve significantly longer operating distances while maintaining the superior chromatic aberration performance of apochromatic lenses.



*Brightfield observation (EPI) objective lens

Model	Magnification	NA	Working Distance(mm)
TU Plan Apo EPI (brightfield type)	50x	0.8	2.0
	100x	0.9	2.0
	150x	0.9	1.5
TU Plan Apo BD (brightfield/ darkfield type)	50x	0.8	2.0
	100x	0.9	2.0
	150x	0.9	1.5

Other lenses

Lenses with correction mechanism

CFI L Plan EPI CR Series

EPI 20x/50x/100x




Model	Magnification	NA	Working Distance(mm)	Glass Thickness Correction Range (mm)
CFI L Plan EPI CR	20x	0.45	10.9-10.0	0-1.2
CFI L Plan EPI CR	50x	0.7	3.9-3.0	0-1.2
CFI L Plan EPI CRA	100x	0.85	1.2-0.85	0-0.7
CFI L Plan EPI CRB	100x	0.85	1.3-0.95	0.6-1.3

Microscope camera


Digital Sight 1000

Equipped with a 2 megapixel CMOS image sensor, it can capture full HD microscope images. By connecting a microscope to this camera and HDMI monitor, movies and images can be captured and saved onto a pre-inserted SD card in the camera.

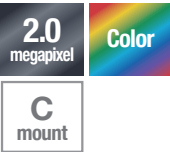


Digital Sight 100

Combined with industrial microscopes, the camera delivers 6.5-megapixel resolution (2944x2208 pixels). HDMI monitor output enables on-site observation without a PC.

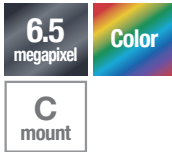

Digital Sight 10

This high-resolution camera captures both color and monochromatic images at up to 6,000 x 3,984 pixels. This enables the wide range of images to be captured and then many of them to be stitched together making a single and large combined image.

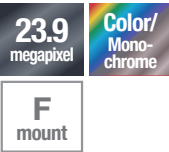

Max Frame Rate	30 fps (1920x1080)	60 fps (1600x900)	55 fps (2000x1328)
Max Recordable Pixels	1920x1080	2944x2208	6000x3984


*Digital Sight 100, standalone, delivers up to 17.7-megapixel resolution (4864 x 3648 pixels).

Imaging software NIS-Elements

Using a desktop PC





Using a desktop PC / tablet PC

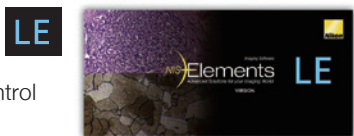

Free software that allows intuitive control of microscope cameras from the PC. Supports Wi-Fi connectivity when used with the Digital Sight 100.


Image Stitching

Stitches together images acquired from multiple fields of view to create one image.

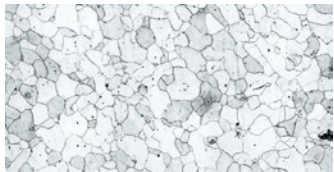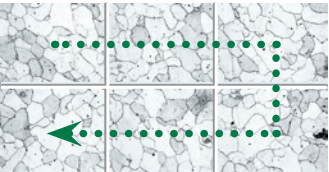



EDF (Extended Depth of Focus)

Create a single, all-in-focus image from images of differing focus.

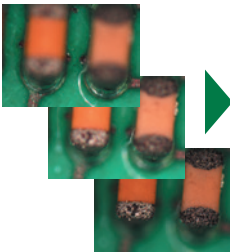

A wide variety of tools

Enables the conducting of simple measurements on images, with input of lines and comments. These can also be written onto and saved with the image, and measurement data can be output.

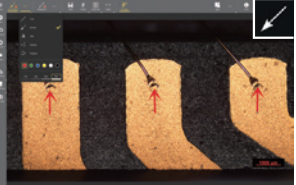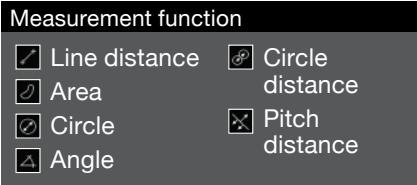

Scene Mode

Ten camera setting patterns for optimal color reproduction and contrast for each microscope light source, observation method and type of sample, as well as custom settings, can be selected.

• Wafer/IC	• Metal, Ceramic/Plastic
• Circuit board	• Flat Panel Display

Information detection and control of accessories are not available with NIS-Elements LE, L, and F. Please use NIS-Elements D for these functions.

Wafer loader NWL200

Combined with the NWL200 wafer loader, the ECLIPSE L200N meets requirements for wafer inspections.

Support for ultra-thin 100 μm wafers

- NWL200 series provides levels of safety and reliability that meet all requirements for inspection of the latest wafers.


Improved operability and high throughput

- Setting conditions, such as sampling and inspection patterns, and checking the operating status and content of errors can easily be done with the large LCD panel
- Comprehensive file management functions for carriers and samples are useful for automating inspections
- Exceptionally fast elevator, and the loading and unloading of wafers with complete precision by the multi-arm system all contribute to an efficient wafer transfer and exchange




Dimensional diagram (Unit: mm)

ECLIPSE L300N/L300ND

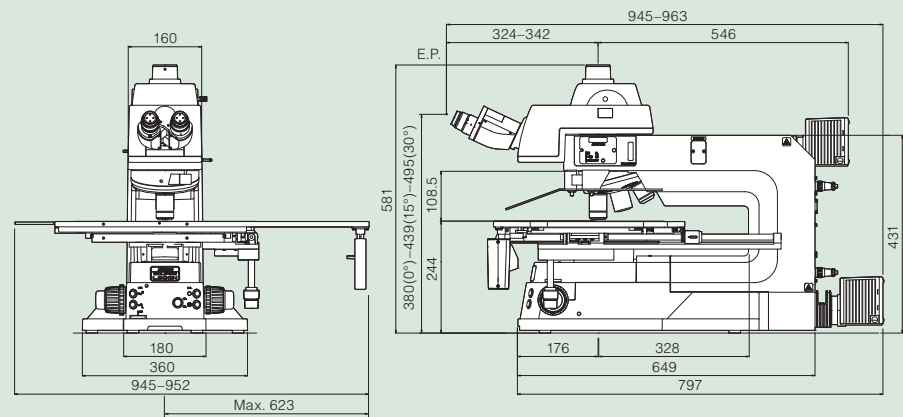


ECLIPSE L200N/L200ND

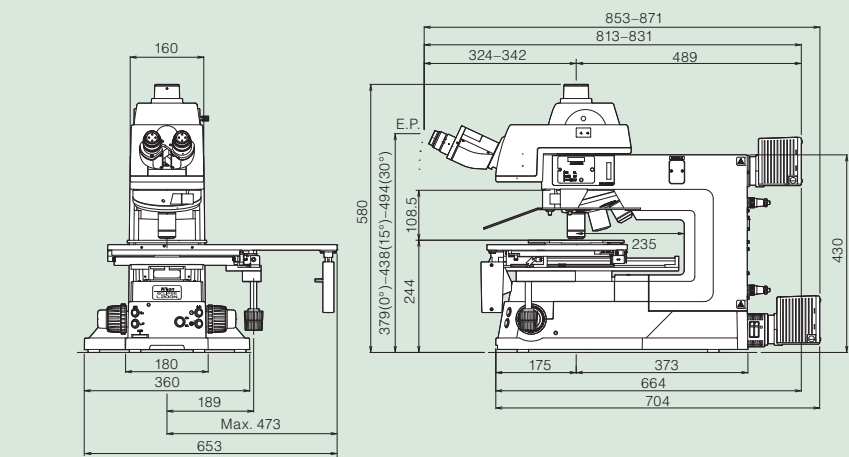


System diagram

ECLIPSE L300N/L300ND

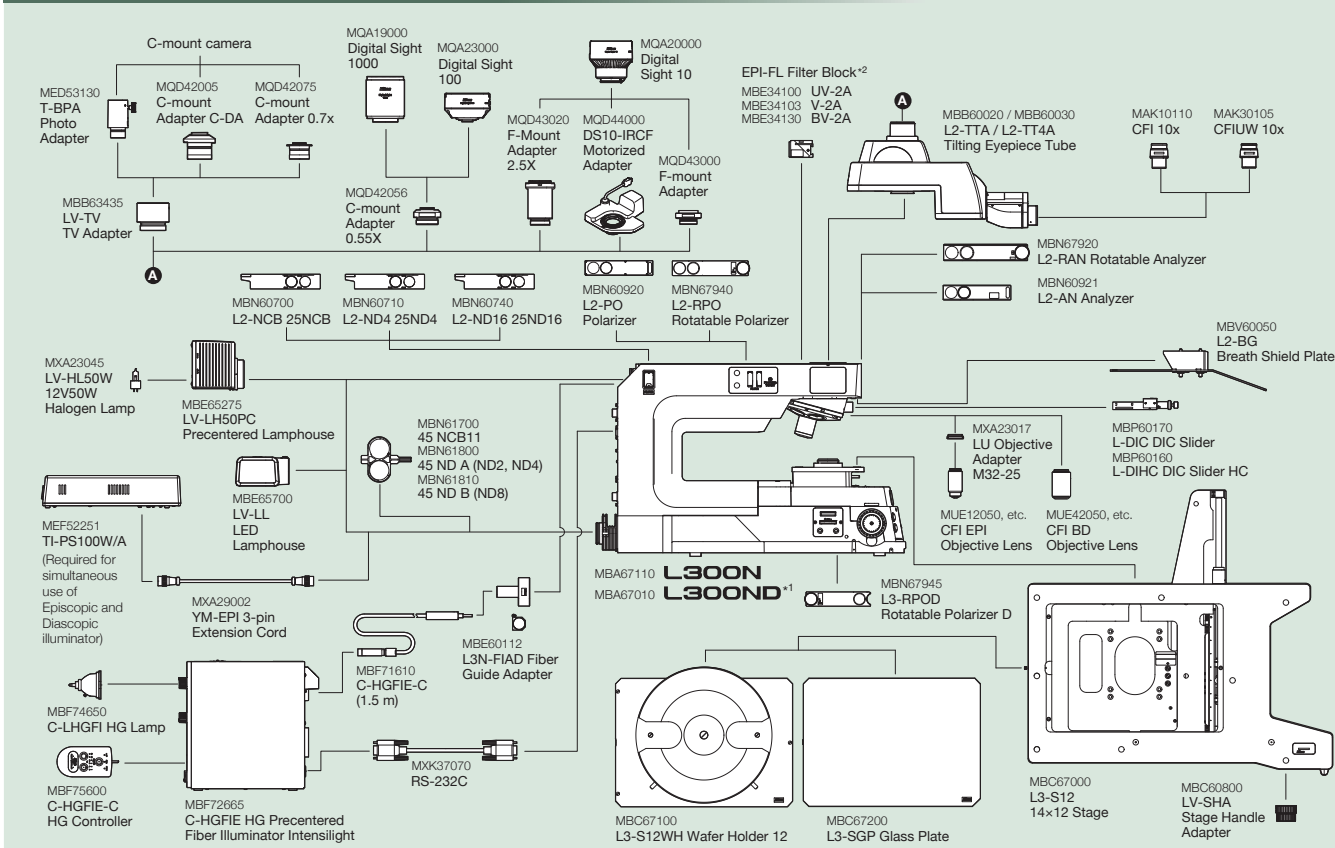

ECLIPSE L200N/L200ND

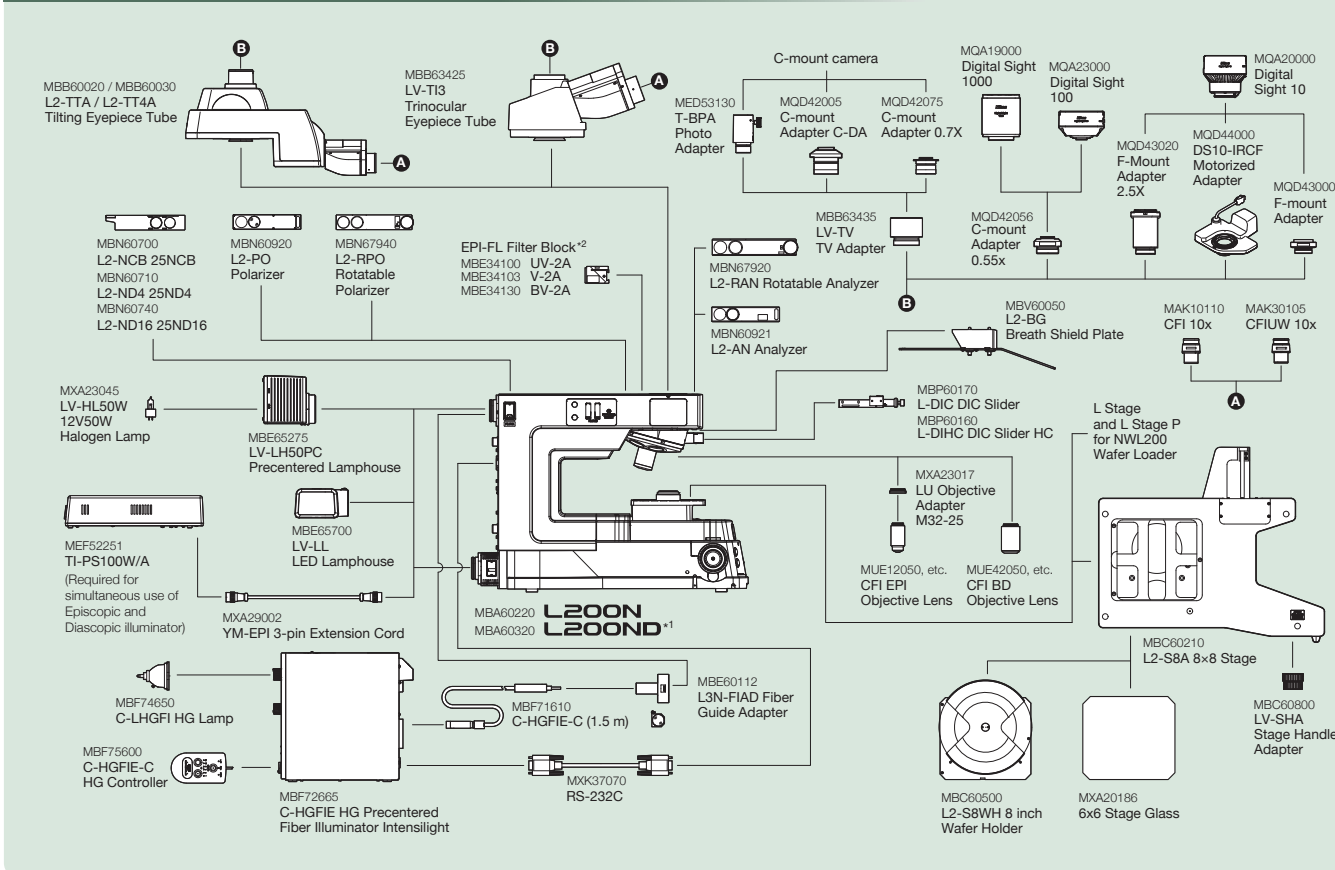

*1 Diascopic illumination available only for L300ND and L200ND
*2 Epi-fluorescence observation available only for L300ND/L300N/L200ND

Specifications

		ECLIPSE L300N	ECLIPSE L200N	ECLIPSE L300ND	ECLIPSE L200ND
Illumination type		Episcopic		Episcopic/Diascopic	
Main body		Power sources for motorized control built in Motorized control for nosepiece, Light intensity control, Aperture diaphragm control			
Nosepiece		Motorized universal sextuple nosepiece			
	Centering Function	Yes	—	Yes	—
	EPI/DIA changeover	—	—	Yes	
Focusing mechanism	Cross travel	29 mm			
	Coarse	12.7 mm per rotation (torque adjustable, refocusing mechanism provided)			
	Fine	0.1 mm per rotation (in 1 μm increments)			
Episcopic illuminator		12V-50W halogen lamp light source built in, LV-LL LED Lamphouse Motorized aperture diaphragm (centerable), Fixed field diaphragm (with focus target) Pinhole slider (optional), Four ø25 mm filters (NCB11, ND16, ND4), Polarizer and Analyzer can be mounted Observation methods: Brightfield, Darkfield, Simple polarizing, DIC, Epi-fluorescence* (*L300N/L300ND/L200ND only)			
Diascopic illuminator		—		12V-50W halogen lamp light source built in, LV-LL LED Lamphouse Aperture diaphragm built in LWD condenser built in	
Interface		USB x 1, RS232C (for Intensilight) x 1			
Eyepiece tubes		Tilting Trinocular Tube L2-TTA (tilt angle: 0-30°) FOV: 22/25; Beam split ratio 100:0/0:100 Tilting Trinocular Tube L2-TT4A (tilt angle: 0-30°) FOV: 22/25; Beam split ratio 100:0/20:80 Trinocular Tube LV-TI3 (Erected image) FOV: 22/25; Beam split ratio 100:0/0:100			
Eyepieces		CFI eyepiece lens series			
Objective lenses		CFI60-2/CFI60 system			
Stages		L3-S12 14 x 12 stage	L2-S8A 8 x 8 stage	L3-S12 14 x 12 stage	L2-S8A 8 x 8 stage
	Stroke	354 x 302 mm	205 x 205 mm	354 x 302 mm	205 x 205 mm
	Diascopic observation range	354 x 268 mm	150 x 150 mm	354 x 268 mm	150 x 150 mm
		Coarse/Fine-movement changeover possible, Fixed-position X-Y fine-movement controls, LV-SHA stage handle adapter attachable			
Antistatic mechanism		1000-10 V, within 0.2 sec			
Power consumption		1.2 A/90 W			
Weight (approx.)	Body only	38 kg	31 kg	40 kg	31 kg
	With L2-S8A 8 x 8 stage and L2-TTA eyepiece tube	—	45 kg	—	44 kg
	With L3-S12 14 x 12 stage and L2-TTA eyepiece tube	57 kg	—	59 kg	—

Specifications and equipment are subject to change without any notice or obligation on the part of the manufacturer. June 2026 ©2010-2026 NIKON CORPORATION

N.B. Export of the products* in this brochure is controlled under the Japanese Foreign Exchange and Foreign Trade Law. Appropriate export procedures shall be required in case of export from Japan.

*Products: Hardware and its technical information (including software)



WARNING

TO ENSURE CORRECT USAGE, READ THE CORRESPONDING MANUALS CAREFULLY BEFORE USING THE EQUIPMENT.



NIKON CORPORATION

1-5-20, Nishioi, Shinagawa-ku, Tokyo 140-8601, Japan

Tel: +81 3 6743 5742

<https://industry.nikon.com/>

ISO 14001 Certified
for NIKON CORPORATION

NIKON METROLOGY EUROPE NV

Interleuvenlaan 86, 3001 Leuven, Belgium

Tel: +32 16 74 01 01

Sales.NM@nikon.com

NIKON METROLOGY UK LTD.

UNITED KINGDOM Tel: +44 1332 811 349

Sales.UK.NM@nikon.com

NIKON METROLOGY SARL

FRANCE Tel: +33 1-60 86 09 76

Sales.France.NM@nikon.com

NIKON METROLOGY GMBH

GERMANY Tel: +49 211 4544 6951

Sales.Germany.NM@nikon.com

NIKON METROLOGY, LLC

12701 Grand River Road, Brighton, MI 48116 U.S.A.

Tel: +1 810 220 4360

Sales.NM-US@nikon.com

NIKON METROLOGY - MEXICO

Sales.NM-US@nikon.com

NIKON PRECISION (SHANGHAI) CO., LTD.

CHINA Tel: +86 21 6841 2050 (Shanghai)

CHINA Tel: +86 10 5831 2028 (Beijing branch)

CHINA Tel: +86 20 3882 0551 (Guangzhou branch)

Web.Nis@nikon.com

NIKON INSTRUMENTS KOREA CO., LTD.

KOREA Tel: +82 2 6288 1900

NIK.Sales@nikon.com

NIKON SINGAPORE PTE LTD.

SINGAPORE Tel: +65 6559 3651

NSG.Industrial-sales@nikon.com

PT. NIKON INDONESIA

INDONESIA Tel: +62 213 873 5005

PTN.Instruments@nikon.com

NIKON SALES (THAILAND) CO., LTD.

THAILAND Tel: +66 2633 5100

NST.Inst@nikon.com

IM-L2606EN

En